




www.splint.co.kr

EZ SPLINT

that can be custom molded by moisture-curing

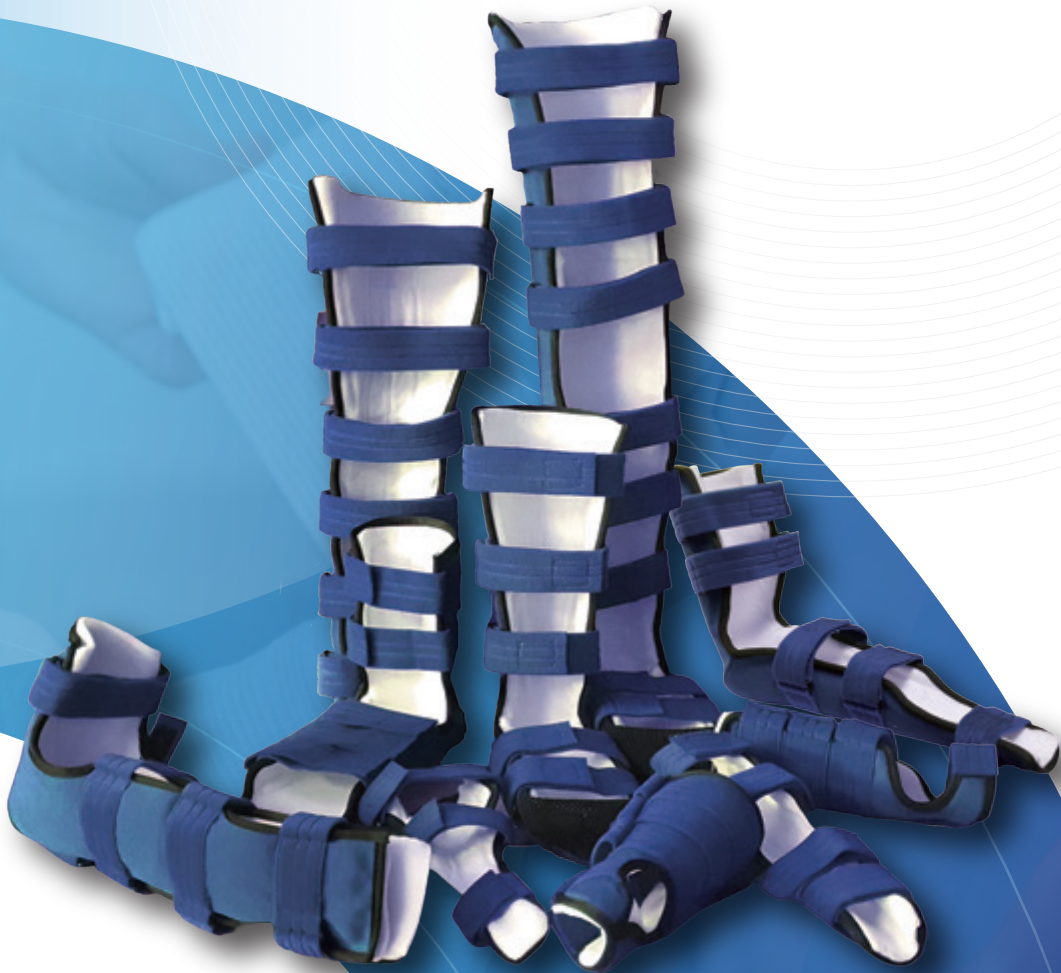

LONG ARM



ELA V110	extra large
ELA V105	large
ELA V100	medium
ELA V98	small plus
ELA V95	small
ELA V90	extra small



SUGAR TONG



ESA VS110	extra large
ESA VS105	large
ESA VS100	medium
ESA VS95	small
ESA VS90	extra small



SHORT ARM

Type 1



ESA V1110	extra large
ESA V1105	large
ESA V1100	medium
ESA V1095	small
ESA V1090	extra small



Type 2



ESA V2110	extra large
ESA V2105	large
ESA V2100	medium
ESA V2098	small plus
ESA V2095	small
ESA V2090	extra small



COCK UP R L

Type 3




ESA V3110	extra large
ESA V3105	large
ESA V3100	medium
ESA V3095	small
ESA V3090	extra small



ULNAR GUTTER

VOLAR

Type 4

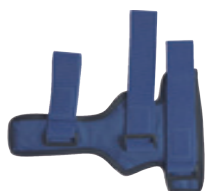

ESA V4110	extra large
ESA V4105	large
ESA V4100	medium
ESA V4095	small
ESA V4090	extra small



THUMB SPICA

LONG LEG



ELL V300	extra large
ELL V290	large
ELL V280	medium plus
ELL V270	medium
ELL V260	small plus
ELL V250	small
ELL V240	extra small



CYLINDER



ELL VC110	extra large
ELL VC105	large
ELL VC100	medium
ELL VC98	small plus
ELL VC95	small
ELL VC90	extra small



SHORT LEG



ESL V300	extra large
ESL V290	large
ESL V280	medium plus
ESL V270	medium
ESL V260	small plus
ESL V250	small
ESL V240	extra small



HALLUX VALGUS



ESL VM110	extra large
ESL VM105	large
ESL VM100	midium
ESL VM95	small
ESL VM90	extra small



Key Features



Anatomically Pre-shaped Design

Conforms naturally to body contours without the need for cutting, heating, or trimming.



Quick & Efficient Application

Significantly reduces application time, improving clinical workflow and treatment outcomes.



Lightweight Yet Rigid

Offers excellent rigidity while maintaining superior comfort and wearability.



X-ray Compatible

Provides excellent radiolucency for easy imaging without removal.



Secure Fit with Velcro Straps

Eliminates the need for additional elastic bandages or external fasteners.



Removable Splint

Can be attached and detached by the patient.


Instructions for Use

1



Open the sealed package.

2

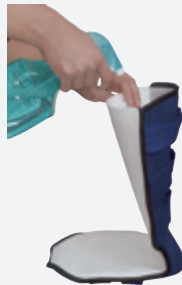

Spread both sides and spray sufficient water.

3



Fasten the Velcro strap through the hook.

4



Mold the splint to fit the body shape. (Takes about 5 minutes)

5



Cover with the foot protector.

6



The splint is now ready for use.



Main Office

2 Jangheung 3-gil, Yangsan-si, Gyeongsangnam-do, Republic of Korea T. +82-55-363-3730 F. +82-55-363-3732

Seoul Office

118, Seongsui-ro, Seongdong-gu, Seoul, Republic of Korea T. +82-2-588-0301 F. +82-55-363-3732

E-mail. cnssplint@gmail.com / cnskorea17@naver.com