

Unique ergonomic design

EZ cast hardened by moisture

DONGHAE MEDICAL
Short Leg



Long Leg



Long Leg(Cylinder)



Long Arm



Short Arm

EZ cast

DONGHAE
MEDICAL

- ✓ Since the connecting part is made of soft plastic, it can be cut with a knife.



- ✓ Excellent fixability since there are supports at the front and rear.

- ✓ Breathability and itchiness relief through the exposed area.

- ✓ Fixed with Velcro so that it can be adjusted to suit the patient's body type.

I Type

Arm

Short Arm



Long Arm



Leg

Short Leg



Long Leg



Long Leg(Cylinder)



I Features **1st** World's first innovative product



Convenience of removal

Since the connecting part is made of soft plastic, it can be cut with a knife, scissor, cast saw. (Patent registered)



Breathability & Cleanliness

It relieves stuffiness, has excellent breathability, and relieves itchiness through the exposed area.



Saves medical procedure time

It is easy to perform a medical procedure since there is no need to bend, twist, or pull as it is made to fit the shape of a leg or an arm.



Excellent strength

It is high-strength and unbreakable since it is made of polyester fabric.



Excellent fixability

It has excellent fixability that allows superior treatment effect since there are supports at the front and rear.



Light weight



X-ray penetration

X-ray can be taken while it is worn.

I Precautions


- Store in a dry and cool area at room temperature (about 24°C/ 75°F)
- The curing time can be adjusted by the time soaked in water and the temperature of the water
- Disposable product
- Shelf life : 3 years

I Size Guide

- L Large** Adults over 90kg (200lb)
- M Medium** General adults
- S Small** Adult female, Adolescents (Ages 13 to 18)
- ES Extra Small** Children (under 12 years old)

Type	Model				Unit
Short Arm	EZSA 1105 L	EZSA 1100 M	EZSA 1095 S	EZSA 1090 ES	10/Box
Long Arm	EZLA 105 L	EZLA 100 M	EZLA 95 S	EZLA 90 ES	5/Box
Short Leg	EZSL 105 L	EZSL 100 M	EZSL 95 S	EZSL 90 ES	5/Box
Long Leg	EZLL 105 L	EZLL 100 M	EZLL 95 S	EZLL 90 ES	5/Box
Long Leg(Cylinder)	EZLL C105 L	EZLL C100 M	EZLL C95 S	EZLL C90 ES	5/Box

I Certificate

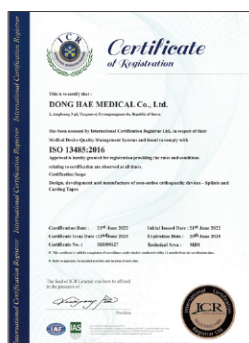

ISO 13485:2016



USA Patent



Japan Patent



China Patent



Korea Patent

Product Manual

1



Wear Stokinet.

2



Spray a sufficient amount of water on the product.

3



Connect the Velcro after wearing the product.

4



Spray a sufficient amount of water.

5



Rub and mold the product while wearing gloves.

6



Completed Cast.

How To Remove



▶ Cut the connecting part using a knife, scissor, cast saw.



EZ CAST

You can watch a video on the procedure and removal method of EZ CAST



Main Office

2 Jangheung 3-gil, Yangsan-si, Gyeongsangnam-do, Republic of Korea T. +82-55-363-3730 F. +82-55-363-3732

Seoul Office

118, Seongsui-ro, Seongdong-gu, Seoul, Republic of Korea T. +82-2-588-0301 F. +82-55-363-3732

E-mail. cnssplint@gmail.com / cnskorea17@naver.com